
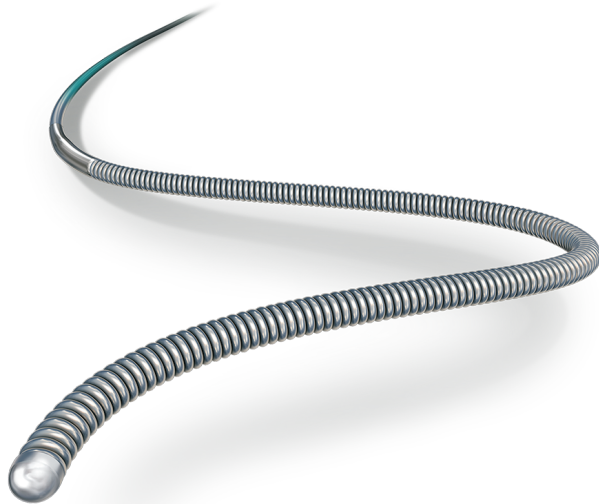


**VASMED TECHNOLOGIES LTD. FZC**  
Facility No. 116, Fujairah, Free Zone  
Fujairah, United Arab Emirates  
Tel.: +971 (0) 9 222 88 62  
Fax: +971 (0) 9 222 88 63  
Email: [sales@vasmed.ae](mailto:sales@vasmed.ae)  
[www.vasmed.ae](http://www.vasmed.ae)

Avery label J8162, 63 (9.91 x 3.81 cm)

**AcceX**  
0.014" PTCA Guidewire

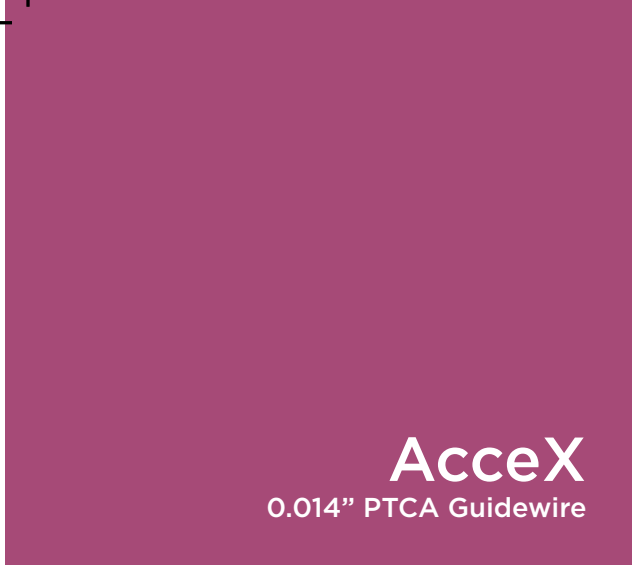


The Vasmed PTCA Guide wire Series are intended to facilitate the placement of devices during diagnostic and interventional procedures. The guide wires have a distal section made of hydrophilic coated coil, the proximal section is made with a SS hardened core wire. The core wire is fixed at the proximal end and extends to tip of the distal end. The proximal section of the core wire is PTFE coated. The distal segment of the core wire has a profiled taper. The core wire enhances the stiffness of the guide wire. A coil is provided at the distal section of the core wire up-to 30 cm, which shall provide enough flexibility to the distal section of the guidewire. The tip of the guidewire is constructed with a Platinum alloy marker for visualization under fluoroscopy.

The Vasmed PTCA guide wire series are offered in 0.014 inch outer diameter with different tip configuration (straight or J tip) with various support of Floppy, Normal and extra support and are available in length 185cm.

INTENDED USE OF DEVICE

The PTCA Guide Wire is used to facilitate the placement of devices during diagnostic and interventional procedures. The guide wires are intended to facilitate catheter placement and exchange during diagnostic or interventional procedures, where increased support, distal flexibility and low surface friction of the guide wire is needed.

KEY PARAMETERS
Unicore design
Multiple core support options
Hydrophilic Distal shaft
Excellent Maneuverability
Good Support for device delivery
Excellent trackability
Good Workhorse wire

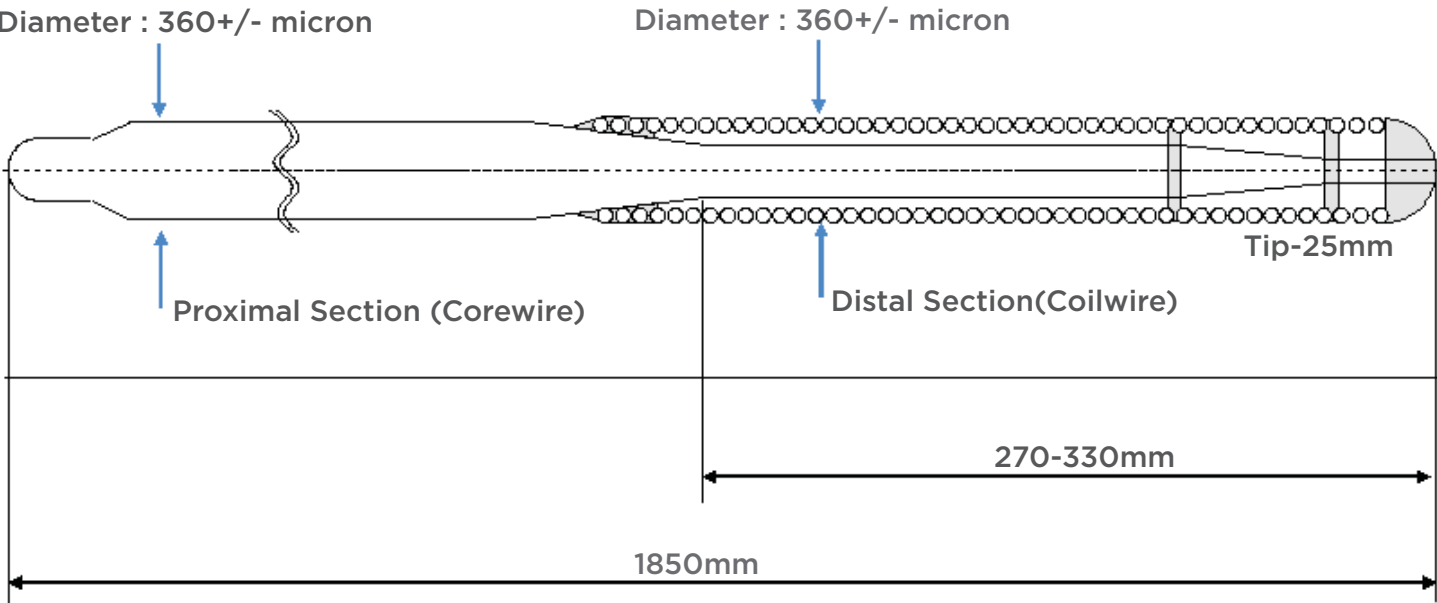

Proximal Shaft	Stainless Steel
Distal Shaft	Stainless Steel
Tip	Platinum Alloy
Coating - Proximal	PTFE coating
Coating - Distal	Hydrophilic coating

Sl. No	Part No.	Wire Dia.	Wire length	Tip profile	Tip length.	Support Type.
1	60000008	0.014"	1850	Straight	25mm	Floppy
2	60000009	0.014"	1850	J	25mm	Floppy
3	60000010	0.014"	1850	Straight	25mm	Normal
4	60000011	0.014"	1850	J	25mm	Normal
5	60000012	0.014"	1850	Straight	25mm	Extra Support
6	60000013	0.014"	1850	J	25mm	Extra Support